



The Grid 2 – creating word lists for prediction

In addition to learning new words for word prediction as you type, The Grid 2 also gives you the option of adding additional words to the prediction list in other ways. This instruction sheet will explain how to add a list of words to be predicted. This gives you the advantage of prioritising certain words for prediction.


1. The first step is to type the words you would like to add into a word document. One word per line. Then save the document as a text or .txt document. To do this select Save as and choose .txt or plain text file from the list of file types.
2. Enter editing mode (**Edit→Editing mode**) and select **Edit Grid Settings**.

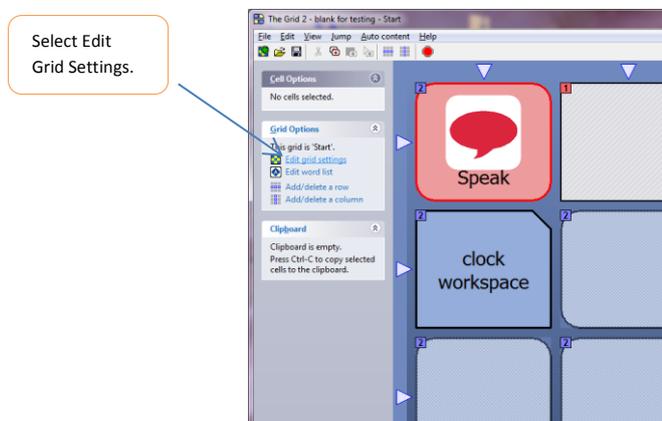

3. Select the **Auto content** tab. Choose the most appropriate option for your needs. Most commonly this will be **Show words from the word list first, then use the full dictionary**. This will ensure that words in your word list are predicted first and then other words will be predicted in subsequent cells.
4. Select **Edit Word List**

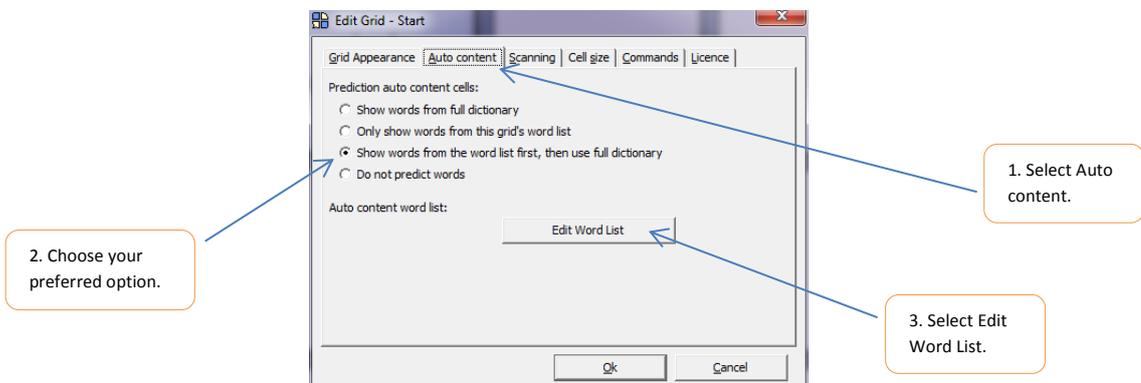

5. Select **File→ Add words from a text document**. Then choose the text document where you have saved the words you want to add. Choose whether you would like the Grid 2 to find pictures to match the words. Select **Save and Close**.

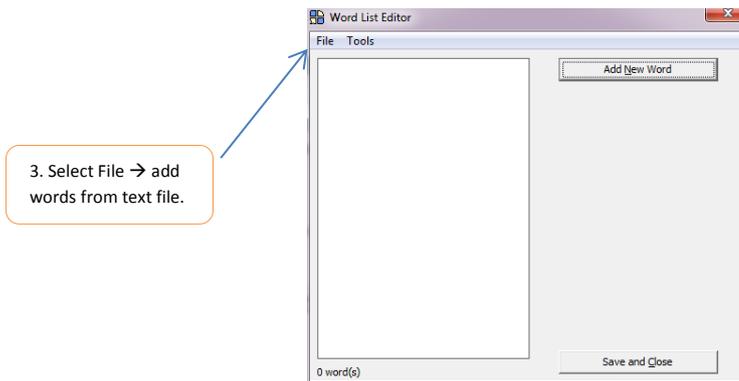

6. If your page doesn't already have prediction cells you can include them. You can select more than 1 cell at a time by dragging the mouse around the cells you want to select. Select **Auto content cells** from the left hand side and select **prediction**.



7. Exit editing mode by selecting the **stop** sign at the top of the screen. Don't forget to save your changes.

Your new words will now be predicted on this page. You can repeat this process with the same and/or different lists of words on additional pages in your grid set.